


Fixed Height Hospital Beds

LABAIRE



Four, Two & One Section Beds

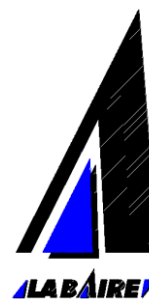

Operation

The 2100 series Labaire Fixed Height Hospital Beds are designed for patient comfort and nursing assistance. The fixed height limits the positioning however the tilt and section adjustment allows the patient and nursing staff a practical means of adjusting the bed to suite their requirements. The options without tilt progressing to a single section progressively limits adjustment of the bed.

Applications

The 1900 series beds are designed for general and specialist applications, which include medical and surgical nursing, orthopaedics and paediatrics.

Design Features 2 & 4 x Section Bed

The gas assisted back rest maybe set at any suitable angle and released. The two way tilt allows a 12 degree foot-up and 6 degree foot-down. The four points of support provide complete stability at any tilt degree and patient positioning. The bed is mounted on two free and two brake 125mm castors. The bed is engineered to withstand robust handling and, with its superior baked paint finish and buffer wheels, will last for years. The head and foot ends can be custom designed with the Hospital Group logo incorporated.

Design Features 2 & 4 x Section Bed

The gas assisted back rest maybe set at any suitable angle and released. The bed is mounted on two free and two brake 125mm. The bed is engineered to withstand robust handling and, with its superior baked paint finish and buffer wheels, will last for years. The head and foot ends can be custom designed with the Hospital Group logo incorporated.

Design Features 1 x Section Bed

The bed is mounted on two free and two brake 125mm. The bed is engineered to withstand robust handling and, with its superior baked paint finish and buffer wheels, will last for years. The head and foot ends can be custom designed with the Hospital Group logo incorporated.

Accessories

- Safety sides, folding or free standing
- Bed fitted IV pole with adjustable height
- Lifting Pole
- Bed stripper
- Lithotomy Stirrups and leg rests
- Bottle and cylinder holders
- Orthopaedic attachments

Technical Data

Length	(mm)	2220	Paediatric	1675
Width	(mm)	910		
Height	(mm)	440 to 725		
Mass	(kg)	98		
Working load	(kg)	200		



See our separate hospital bed accessory brochure for details on all accessories available.



Labaire

sales@labaire.co.za

Cell: 082 211 5118

Distributed by: